

Evacuation Chair



The Stryker Evacuation Chair is essential for the evacuation of disabled or injured persons weighing up to 500 pounds, from multilevel facilities in emergency situations. The patented Stair-TREAD tracks allow easy evacuation down stairs by a single operator, for passengers weighing up to 200 pounds. The extendable handle allows the operator to have a forward-facing, upright posture for balance, line-of-sight and ergonomic safety. Foot end lift handles also extend, allowing two people to lift an evacuee over obstacles or up stairs.

Evacuation Equipment

- **Foldaway foot support provides additional passenger security and stability**
- **Contoured seat for comfort**
- **Foot end lift handles extend for ergonomic lifting**
- **Head end lift handles extend for upright lifting posture**
- **Folds and stores in a small space**
- **The patented Stair-TREAD™ system controls descent in proportion to the passenger's weight**
- **Wheels provide maneuverability in narrow passageways**

500^{lb}
capacity

15 YEAR
Warranty



Evacuation Chair

Features and Benefits

- **Innovative Stair-TREAD system** allows easily controlled stair descent
- **Easy set-up and instruction pictograms** on seatback
- **Color-coded controls** for faster operation
- **Extendable head end and foot end lift handles** for ergonomic lifting
- **Extendable upper-control handle** for upright operator posture and unobstructed line of sight
- **Positive action locking mechanism** for passenger security
- **Dual wheel locks** for more secure transfers
- **Oversized rear wheels with sealed bearings** for smoother transports
- **Compact storage size conserves space** in corridors and closets
- **Lightweight, rugged aluminum construction** for low effort transports
- **Power washable** for easy cleaning
- **Grease-free maintenance** for ease of care
- **Restraints¹ with Fastex buckles** for passenger security
- **Molded hand grips fit hands** for better ergonomics
- **Head support allows evacuation of passengers with compromised head control**
- **Two-piece molded ABS seat** for greater passenger comfort

Specifications²

Evacuation Chair

Model Number	6254
Height	45 in. (114.3 cm)
Width	20.5 in. (52 cm)
Depth³	28 in. (71 cm)
Folded Depth	11 in. (28 cm)
Weight	34 lb. (15 kg)
Maximum Load^{4,5}	500 lb. (227 kg)

Accessories




Wall Mount
6254-050-000



Chair Cover
6254-055-000



Foot Support
6254-029-000

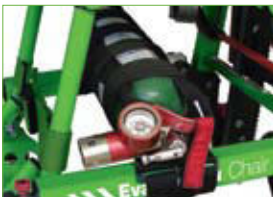

O₂ Holder
6250-140-000



Storage Cabinet
6254-002-000



Alarm
6254-045-000

Evacuation Equipment

3800 E. Centre Ave.
Portage, MI 49002 U.S.A.
t: 269 329 2100 f: 269 329 2315
toll free: 800 784 4336

- ¹ Passenger restraint set includes lap, chest and ankle restraints.
- ² Dimensions are measured from the outermost edges of the main frame. Specifications are rounded to the nearest whole number. Conversions are calculated before rounding.
- ³ Depth dimensions are measured with extendable handles retracted.
- ⁴ Operators must consider the weight of the passenger, equipment and accessories when determining the total weight on the product.
- ⁵ Maximum weight capacity is total weight distributed in accordance with basic human anatomy.
- Stryker reserves the right to change specifications without notice.

Optional



In-service Video
6253-001-167



Aggressive Stair-TREAD Belt[†]
Available as an option

*The Aggressive Stair-TREAD Belt is designed to improve the grip of the Stair-TREAD system on steps that have a smooth surface with a rounded edge. To determine if the Aggressive Stair-TREAD Belt is appropriate for your application, please refer to the Aggressive Stair-TREAD Belt Specification Sheet at www.strykerevacuation.com or contact your local Evacuation Consultant.